



Get Clarity On Generics

Cost-Effective CT & MRI Contrast Agents



FRESENIUS
KABI

WATCH VIDEO

AJNR

Patient Radiation Dose Management in the Follow-Up of Potential Skin Injuries in Neuroradiology

E. Vano, J.M. Fernandez, R.M. Sanchez, D. Martinez, L. Lopez Ibor, A. Gil and C. Serna-Candel

This information is current as of August 28, 2025.

AJNR Am J Neuroradiol 2013, 34 (2) 277-282

doi: <https://doi.org/10.3174/ajnr.A3211>

<http://www.ajnr.org/content/34/2/277>

Patient Radiation Dose Management in the Follow-Up of Potential Skin Injuries in Neuroradiology

E. Vano, J.M. Fernandez, R.M. Sanchez, D. Martinez, L. Lopez Ibor, A. Gil, and C. Serna-Candel



ABSTRACT

BACKGROUND AND PURPOSE: Radiation exposure from neurointerventional procedures can be substantial, with risk of radiation injuries. We present the results of a follow-up program applied to potential skin injuries in interventional neuroradiology based on North American and European guidelines.

MATERIALS AND METHODS: The following guidelines approved in 2009 by SIR and CIRSE have been used over the last 2 years to identify patients with potential skin injuries requiring clinical follow-up: peak skin dose >3 Gy, air kerma at the patient entrance reference point >5 Gy, kerma area product >500 Gy $\cdot$ cm², or fluoroscopy time >60 minutes.

RESULTS: A total of 708 procedures (325 in 2009 and 383 in 2010) were included in the study. After analyzing each dose report, 19 patients (5.9%) were included in a follow-up program for potential skin injuries in 2009, while in 2010, after introducing several optimizing actions and refining the selection criteria, only 4 patients (1.0%) needed follow-up. Over the last 2 years, only 3 patients required referral to a dermatology service.

CONCLUSIONS: The application of the guidelines to patient radiation dose management helped standardize the selection criteria for including patients in the clinical follow-up program of potential skin radiation injuries. The peak skin dose resulted in the most relevant parameter. The refinement of selection criteria and the introduction of a low-dose protocol in the x-ray system, combined with a training program focused on radiation protection, reduced the number of patients requiring clinical follow-up.

ABBREVIATIONS: CIRSE = Cardiovascular and Interventional Radiology Society of Europe; DAP = dose-area product; DOLIR = dose on-line for interventional radiology; ICRP = International Commission on Radiologic Protection; SIR = Society of Interventional Radiology

In recent years, the benefits and complexity of interventional procedures have increased substantially and the radiation dose to patients and staff is a matter of concern.¹⁻³ Among medical specialties, interventional neuroradiology is of utmost benefit to patients, but it uses high radiation doses in a significant number of procedures. Radiation injuries in patients may occur and should

be listed as a risk in the informed consent.^{2,4} In their quality programs, interventional radiology units include the evaluation of patient radiation doses and the criteria for including patients in a clinical follow-up as part of postprocedural care when relevant dosimetric parameters indicate potential skin radiation injuries.^{2,4}

In 2000, the ICRP published a set of recommendations on how to avoid radiation injuries in medical interventional procedures.² These recommendations have been developed in several guidelines and adopted by interventional radiology medical societies.⁴

The European Directive on Medical Exposures 97/43/Euratom⁵ requires Member States of the European Union (article 9) to use appropriate radiologic equipment submitted to quality assurance programs and to assess patient doses. Preventing high-dose exposures with diagnostic equipment is another requirement (article 11).

Many patients are still neither being counseled on radiation risks nor followed up when radiation doses from difficult procedures may lead to injury. In some interventional procedures, skin doses to patients exceed the threshold for deterministic effects.⁶

Received February 26, 2012; accepted after revision April 29, 2012.

From the Medical Physics Service and Radiology Department (E.M., J.M.F.), Instituto de Investigación Sanitaria Hospital Clínico San Carlos and Complutense University, Madrid, Spain; Medical Physics Service (R.M.S., D.M.) and Interventional Neuroradiology Unit, Diagnostic Radiology Service (L.L.I., A.G., C.S.-C.), Instituto de Investigación Sanitaria Hospital Clínico San Carlos, Madrid, Spain.

The authors acknowledge the support of the Spanish grant SAF2009-10485 (Ministry of Science and Innovation).

Previously presented in part as an electronic poster at: CIRSE Congress, September 10–14, 2011; Munich, Germany.

Please address correspondence to E. Vano, Medical Physics Service, Hospital Clínico San Carlos, P. Martín Lagos St, 28040 Madrid, Spain; e-mail: eliseov@med.ucm.es

Indicates open access to non-subscribers at www.ajnr.org

<http://dx.doi.org/10.3174/ajnr.A3211>

Several evaluations of cumulative skin dose have been made in interventional neuroradiology procedures.⁷⁻¹³ Radiation-induced skin injuries can sometimes occur after a clinically complex procedure, but may also, on other occasions, result from the use of inappropriate equipment or poor operational techniques. The ICRP states that acute radiation doses may cause erythema at 2 Gy and delayed skin necrosis at 12 Gy.² Other more recent reports have analyzed the different dose thresholds for radiation injuries in interventional procedures.⁶

The term “cumulative skin dose,” introduced by the ICRP,² is used in this article as the value of the dose in air, accumulated at the entrance of the patient during the whole interventional procedure. It is obtained from the air kerma accumulated at a specific point in space relative to the fluoroscopic gantry (called “patient entrance reference point” by the International Electrotechnical Commission)¹⁴ during the full procedure, and by taking into account the calibration factor of the ionization transmission chamber of the x-ray system and the attenuation of the table and head support and the increase caused by the backscatter factor.

The accumulated air kerma displayed by the x-ray system and presented in the patient dose reports includes neither the backscatter nor the attenuation of the table, nor the patient’s head support. The dose values reported by the x-ray systems also need to be corrected by the corresponding calibration factor of the internal transmission ionization chamber. The “cumulative skin dose” is not the same quantity as the “peak skin dose” (ie, the highest dose at any portion of a patient’s skin during a procedure).

The current standard for the x-ray systems used in interventional radiology¹⁴ requires supplying information on the kerma area product (also known as DAP; usually measured in Gy · cm²) and cumulative air kerma (or cumulative dose) at the patient entrance reference point (usually displayed in mGy). The cumulative air kerma is calculated by the x-ray system (free in air) at 15 cm from the isocenter in the focus direction (the “reference point”), and its value, once corrected by the calibration factor of the internal transmission ionization chamber, is similar to the cumulative skin dose ± 10% in most of the systems (considering the increase caused by the backscatter of the patient and the attenuation in the table and mattress or patient head support).^{15,16}

The ICRP has recommended that maximum cumulative absorbed doses to the skin approaching or exceeding 1 Gy (for procedures that may be repeated) or 3 Gy (for any other procedure) should be recorded in the patient record and that there should be follow-up procedures for such cases.²

In 2009, SIR (in North America) and CIRSE published a common “Guidelines for Patient Radiation Dose Management,” previously adopted by both societies.⁴ The document introduced the term “significant radiation dose” as a selected threshold value used to trigger additional dose-management actions. In our hospital, the procedure adopted in the quality assurance program is in accordance with the postprocedural care included in these guidelines, especially concerning dose documentation and patient follow-up.

Several articles have been published on radiation skin injuries in interventional procedures,¹⁷⁻²² but very few, if any, give details on the percentage of patients exceeding the trigger values and requiring clinical follow-up after potential radiation injuries.

The aim of this article is to present the results of the application of the SIR-CIRSE guidelines in interventional neuroradiology and report the percentages of patients requiring clinical follow-up due to potential radiation injuries.

MATERIALS AND METHODS

The San Carlos Hospital Ethical Committee approved this study under the title “Radiological risks in fluoroscopy guided procedures” (code B-09/20). The x-ray system used in this study was a biplane Allura (Philips Medical Systems, Best, the Netherlands). The system was equipped with a flat imaging detector, with a field of view of 48.4 cm (diagonal dimension) for the frontal C-arm and a flat detector of 25 cm (diagonal dimension) for the lateral C-arm.

The frontal C-arm (48 cm) allows for acquiring rotational and CT-mode series in addition to fluoroscopy and DSA acquisitions (with a large variety of protocols for the different clinical procedures and areas examined). The FOVs available for frontal C-arm detectors are 48, 42, 31, 26, 22, 19, and 15 cm. The lateral C-arm detector (25 cm) has 3 FOVs: 25, 20, and 15 cm. The standard configuration of the system has 3 fluoroscopy modes (all at 15 pulses/second) with added filtration of 1 mm aluminum and 0.9, 0.4, and 0.1 mm of copper for low, normal, and high fluoroscopy modes, respectively.

When the procedures are over, a complete patient dose report is produced by the x-ray system and sent by e-mail (intranet) to the medical physics service, where automatic software (DOLIR)²³ processes the report and, after including the calibration/correction factors, introduces the relevant dosimetric parameters into a data base and into a graphic interface.

The x-ray system is submitted to a quality assurance program and patient dose records are archived individually according to national regulations. The medical physics service is in charge of the quality control of the x-ray and imaging system, dose calibrations, dose records, and patient and staff radiation protection. This service also gives support to the interventional neuroradiology unit to optimize some imaging protocols.


As part of a full optimization program, and as suggested by the medical physics service and the neuroradiology unit, the service engineers prepared a low-dose protocol for patients, implemented during 2010. It basically consisted of a reduction of patient entrance dose rate in fluoroscopy (decreasing the number of pulses/second with a slight increase of the dose per pulse) and of patient entrance dose per image during the DSA series. Details of this protocol are included in the Results section and compared with the standard protocol.

Part of the optimization program also included a 20-hour training course focused on radiation protection for fluoroscopy-guided procedures, corresponding to the “second level of radiation protection training” as recommended by the World Health Organization,¹ ICRP,^{2,24} and the European Guidelines on Education and Training in Radiation Protection for Medical Exposures.²⁵

In addition, DOLIR was implemented at the hospital during 2010²³ as an intermediate solution until the DICOM dose-structured reports are available for interventional radiology systems.

Our biplane neuroradiology x-ray system currently has the

Cerebral embolization 2009 (Dose Area Product)

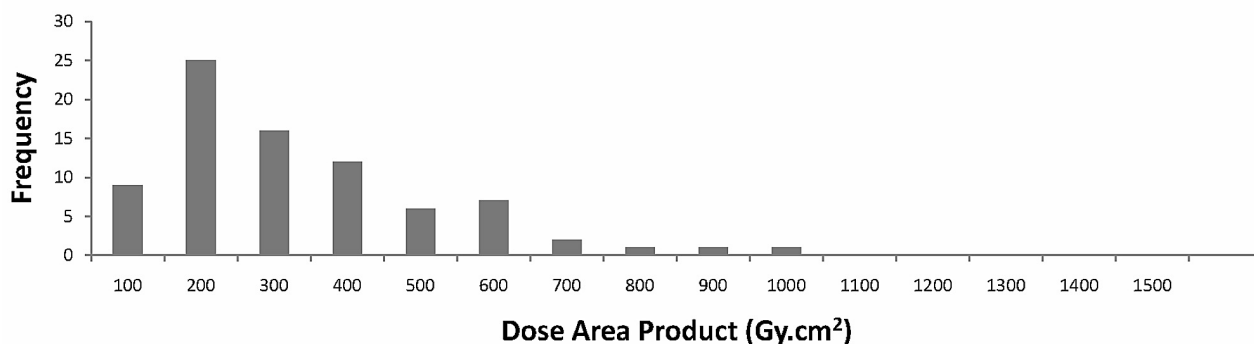

Cerebral embolization 2010 (Dose Area Product)

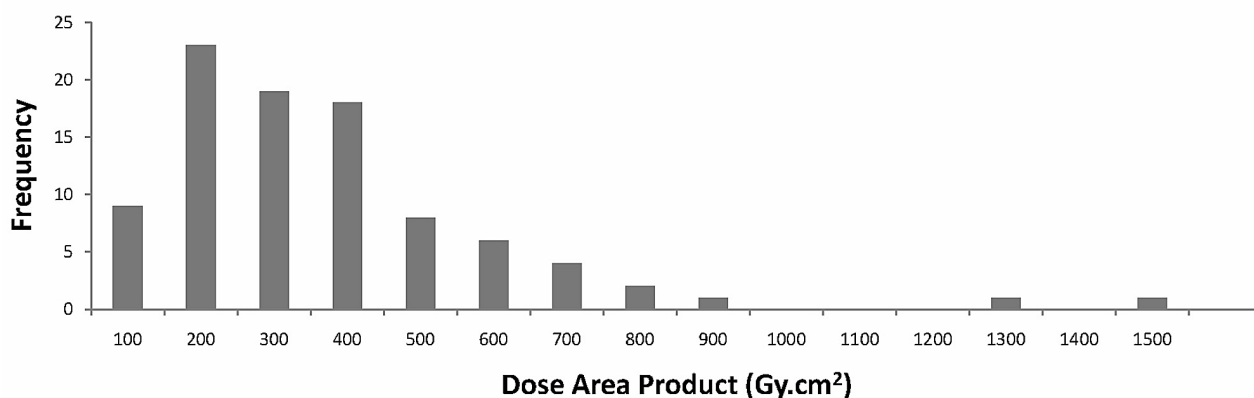


FIG 1. Distributions of dose area product values ($\text{Gy} \cdot \text{cm}^2$) for cerebral embolizations during 2009 and 2010 (frontal and lateral planes together). Median and third quartile values are indicated in Table 1.

capability to export, via e-mail, patient dose reports, including fluoroscopy time, DAP, and cumulative dose at the patient entrance reference point for each of the C-arms at the end of each procedure. Patient dose reports also include details on radiographic techniques, geometry values such as distances and angulations of the C-arms, and the number of frames of all the DSA series and fluoroscopy runs archived.

Dose values are included by the x-ray system as part of the patient dose reports, but in our software, these have been corrected by the corresponding calibration factor (measured as part of the periodic quality controls of the system) to take into account the accuracy of the internal transmission ionization chamber, by the attenuation of the table and the mattress (relevant for the frontal plane) and by the backscatter factor. The values included in the data base are corrected with all these factors.

Trigger levels for a potential patient follow-up were adapted to the values recommended by the SIR-CIRSE guidelines.⁴ These values are peak skin dose >3 Gy, cumulative air kerma at the patient entrance reference point >5 Gy, kerma area product >500 $\text{Gy} \cdot \text{cm}^2$, or fluoroscopy time >60 minutes. The refinement during 2010 consisted of using the cumulative skin dose of each of the planes (frontal and lateral) independently as main trigger levels when one resulted in values >4 Gy, and in taking into account the different angulations of the C-arms during the procedures (included in the patient dose reports) to consider the skin dose distribution.

The medical physics service is alerted when values of DAP or

cumulative doses are over these levels. The dose values used are previously corrected by the calibration of the ionization chamber, the attenuation of the table and head support, and the backscatter factor. The information on patient dose is also displayed in the x-ray room during the procedure.

Alerts are analyzed daily by a senior medical physicist, who decides whether a more detailed analysis of the individual dose reports is appropriate. The information is then transferred to the neuroradiology unit, which, in turn, considers whether it is appropriate to give additional information on potential radiation injuries to the patient or to his or her family and to decide on a clinical follow-up. Part of the analysis consists of searching previous procedures existing in the neuroradiology unit data base and considering whether skin doses are relevant. The neuroradiologist in charge of the patient makes further investigations into the patient's clinical records to check whether he or she could have undergone previous procedures in other hospitals. The clinical follow-up procedure has been approved by the hospital's Quality Assurance and Radiation Safety Committee.

Starting a clinical follow-up, as recommended by the SIR-CIRSE guidelines, means extra work that further adds to the workload of informing patients on the reasons for such a clinical follow-up. Therefore, it is only considered after a careful evaluation of the individual patient dose reports by the medical physics service in the case of dose parameters above trigger levels and especially the peak skin dose for the frontal and lateral planes. These peak skin doses are initially assumed (as a conservative

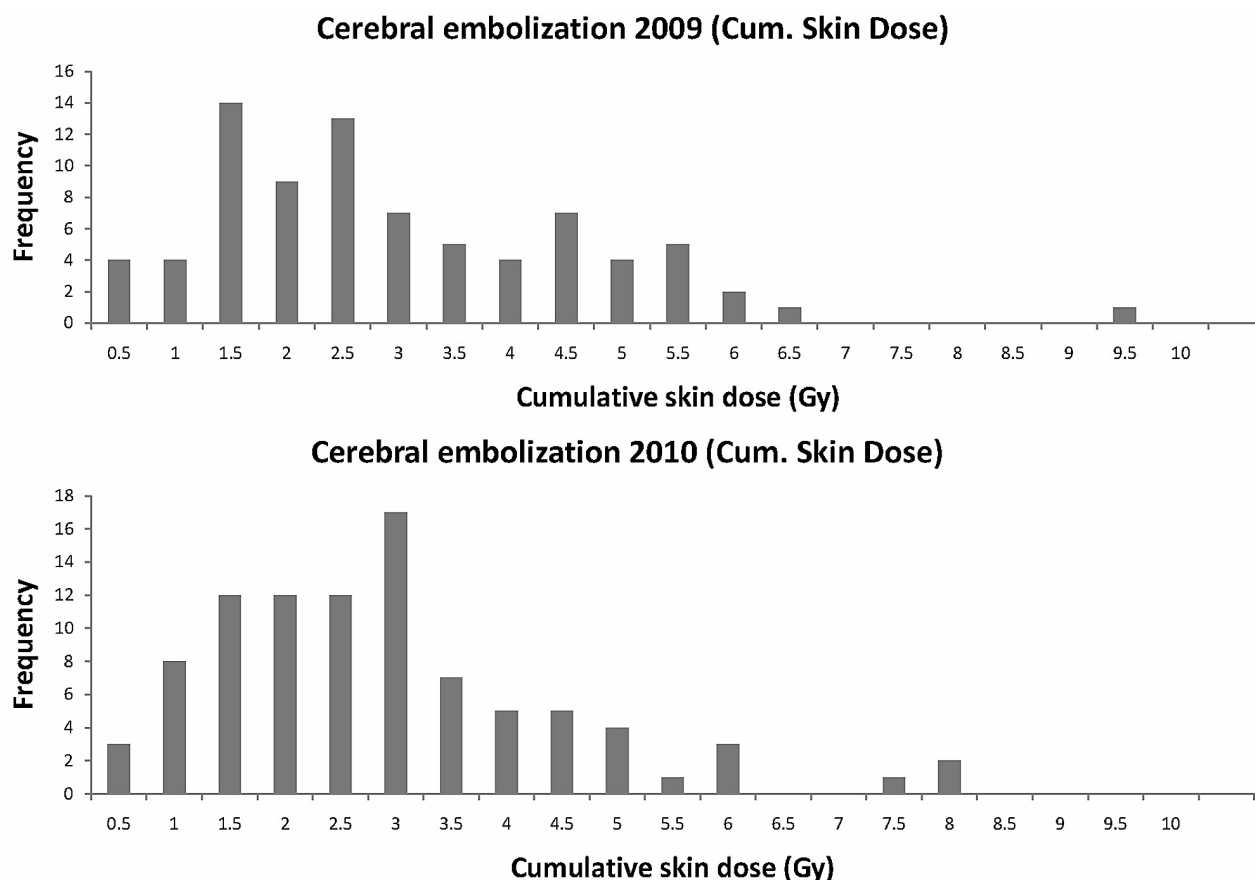


FIG 2. Cumulative skin dose distributions (Gy) for cerebral embolizations during 2009 and 2010 (frontal and lateral planes together). Median and third quartile values are indicated in Table 1.

hypothesis) to be the cumulative skin dose of each of the planes (frontal and lateral), calculated from the cumulative air kerma included in the patient dose report, corrected by the attenuation of the table and head support, and increased by the backscatter factor (for the frontal plane). The different angulations of the C-arms are also considered to discriminate some cases in which the real skin dose distribution could be very different from the above-mentioned initial conservative hypothesis. The final decision concerning whether to admit patients in a clinical follow-up program lies with the neuroradiology unit.

RESULTS AND DISCUSSION

A total of 325 procedures (43% therapeutic, mainly cerebral embolizations) were included in the local data base during 2009, and 383 (40% therapeutic) were included during 2010. From the data base, cases with full patient dose reports were selected for further analysis. It should be noted that, in most interventional suites, patient dose reports are still neither automatically archived nor analyzed. In our center, once the procedure is finished, the operators have to manually send the dose report by e-mail to the medical physics service, where it is processed by the automatic software system DOLIR.²² Some of the dose reports can be lost if the operators fail to send the e-mail. To avoid such mishaps, the manufacturer has been asked to automate e-mails once the procedures are closed.

Distributions of DAP and cumulative skin dose values (frontal

Table 1: Dose values for cerebral embolizations

	DAP (Gy · cm ²)	Cumulative Dose (Gy)	DAP (Gy · cm ²)	Cumulative Dose (Gy)
Year	2009	2009	2010	2010
Sample	80	80	92	92
Median	242	2.4	270	2.5
3rd quartile	386	3.9	392	3.3
Mean	293	2.7	317	2.6
SD	188	1.7	234	1.5

and lateral planes together) for embolization procedures for the years 2009 and 2010 are presented in Figs 1 and 2. During 2009, the trigger level used for clinical follow-up was the total cumulative dose (>5 Gy in frontal and lateral planes together), but it was considered too conservative for cerebral procedures. Therefore, in 2010, the patient was included in the follow-up program when the cumulative skin dose value was >4 Gy (in one of the planes), and the individual analysis of the C-arm angulations indicated that this conservative hypothesis could be adopted (eg, most of the x-ray beam directions irradiate the same skin area).

Table 1 presents the results of the patient dose values for cerebral embolizations during 2009 and 2010. The 2009 median values for DAP and cumulative skin dose were very similar to the 2010 ones. There were no statistical differences for either patient dose distribution ($P = .643$ for DAP, and $P = .988$ for cumulative skin dose with the Mann-Whitney U test). However, third quar-

Table 2: Values of DAP and skin dose for cerebral embolizations

Authors and Reference	Year	Sample	DAP		Skin Dose in			Comments
			(Gy · cm ²)	Mean	SD	Gy	Mean Values ^a	SD
Mamourian et al ⁷	2010	8	NA	NA	NA	0.9	NA	Average skin dose
Moritake et al ⁸	2008	32	NA	NA	NA	1.8	1.3	Peak skin dose
Suzuki et al ⁹	2008	103	257	150	1.1	1.9		Maximum entrance skin dose
D'Ercole et al ¹⁰	2007	42	413	NA	1.2			Maximum skin dose
Moskowitz et al ¹¹	2010	46	NA	NA	10.4	4.2		Cumulative skin dose
Miller et al ¹²	2003	356	320	NA	3.8	NA		Cumulative dose
D'Ercole et al ¹³	2012	82	383	176	NA	NA		No skin dose reported. Median value of DAP: 352
This article	2012	172	305	211	2.7	1.6		Cumulative skin dose. Median value of DAP: 256

Note:—NA indicates not available.

^a See Comments column.

Table 3: Dose reduction in entrance surface dose rate (with backscatter) measured with a PMMA phantom

	Low-Mode Fluoroscopy (mGy/min)	Normal-Mode Fluoroscopy (mGy/min)	DSA (mGy/image)
Normal dose (15 p/s in fluoro.) procedure	9.9	25.2	5.9
Low-dose (7.5 p/s in fluoro.) procedure	7.0	16.6	3.1
Dose reduction % (low dose/normal dose)	30	34	47

Note:—PMMA indicates polymethylmethacrylate. Focus-phantom distance: 67 cm; Focus image detector; 104 cm. 20-cm thickness of PMMA; FOV 31.1 cm. The dose per pulse in the new low-dose fluoroscopy mode has been increased in comparison with the “normal” one to reduce the noise.

tile values were lower in 2010 for cumulative skin dose (3.3 Gy) than in 2009 (3.9 Gy), most likely partly because of the radiation protection training and the application of the low-dose protocol initiated in 2010.

The results on skin doses are difficult to compare with other authors because different metrics have been used. We reported cumulative skin doses, obtained from the cumulative air kerma at the patient entrance reference point and displayed by most of the modern x-ray systems,¹³ while other authors have measured or estimated peak skin doses.^{7–13} Table 2 (adapted by the authors) compares the results reported by several authors. Our cumulative skin doses and DAPs are in the range of other published values.^{7–13} In 2009, 19 patients (5.9% of the total number of procedures included in the data base during that year) were included in a follow-up program for potential skin injuries. That year, 20 procedures (2 in the same patient) resulted in cumulative skin doses (frontal and lateral planes together) higher than 5 Gy. In 2010, after introducing optimization actions (ie, a radiation protection training course and a new “low-dose” protocol in the x-ray system), only 11 procedures resulted in a total cumulative dose higher than 5 Gy.

The low-dose protocol (Table 3) was introduced during 2010 as an alternative to the “standard” protocol. Neuroradiologists select the low-dose protocol (especially in the most complex cases) either for part of the procedure or for the full procedure. A numeric evaluation of image quality with test objects was completed before the clinical use of the low-dose protocol, and the neuroradiologists found the differences (in image quality) with the standard mode satisfactory. Dose reduction in fluoroscopy resulted in 30% savings for the low fluoroscopy mode (the most commonly used) and in 47% for DSA images (measured in the conditions indicated in Table 3).

After implementing the SIR-CIRSE guidelines, and after considering the peak skin dose values independently for the frontal and lateral planes and taking into account the C-arm angulations during the procedures, only 4 patients (1.0% of the total number of procedures included in the data base during that year) were

considered for follow-up. In 2009, 8 patients from the 19 cases (1 of the patients with 2 procedures) included in the follow-up program suffered from temporary alopecia. In 2010, the selection criteria were refined by using the peak skin dose (considered as the cumulative skin dose in each of the frontal and lateral planes and the C-arm angulations), and from the 11 patients with more than 5 Gy (as total cumulative skin dose), only 4 patients were included in the follow-up program (all suffered from temporary alopecia). Only 3 patients (0.4% of the full sample included in the data base) during the 2 years required referral to the dermatology service. In the 2-year period, the maximum patient dose values measured in a single procedure were 1500 Gy · cm² and 9.5 Gy (cumulative skin dose from frontal and lateral planes together).

Part of the optimization program consisted of the routine application of the following dose reduction rules: using the low fluoroscopy mode whenever possible; reducing the fluoroscopy time and the number of DSA images and runs; archiving some of the fluoroscopy runs so as to avoid some DSA series; collimating to the area of interest and using virtual collimation; maintaining the image detector as close as possible to the patient and the x-ray tube as far as possible from the patient's skin (especially relevant for the lateral plane); using magnification only if strictly necessary; and using low-dose protocols, if appropriate.

The median values of these dose distributions are good indexes for drawing comparisons with other similar publications and for moving the project of a future set of “diagnostic reference levels” forward, as recommended by the International Commission on Radiologic Protection²⁶ and the future European Directive on Basic Safety Standards.²⁷

One of the most extensive studies on patient doses for interventional radiology made in the last years is the RAD-IR study in the United States.¹² In this study, the group of neuroradiology interventions for the head included embolization of arteriovenous malformations, aneurysms, and tumors. In the US study, a sample of 382 procedures was collected, with a mean patient dose value of 320 Gy · cm². In our case, mean values for the global 2-year sample (172 procedures identified as “embolizations” in

the data base) resulted in 305 Gy · cm². In any case, median and third quartile values are the best descriptors for highly skewed distributions (Figs 1 and 2).

Limitations of the Study

The clinical complexity of the procedures has not formally been evaluated in this study, but an overview of clinical reports in the neuroradiology unit indicates that the therapeutic procedures were of similar complexity in 2009 and 2010, while the clinical staff of the neuroradiology unit remained unchanged. Despite the calibration factor applied to convert the dose values displayed by the x-ray system to entrance skin dose (including the attenuation of the table and head support and the backscatter factor), inaccuracies can occur in the estimation of the peak skin dose, especially when different C-arm angulations are used during the procedure, when radiation fields overlap, or when tight collimation is used during the procedure causing the real peak skin dose to be lower than the cumulative skin dose in one of the planes.

CONCLUSIONS

The application of the SIR-CIRSE guidelines for patient radiation dose management enabled standardization of the selection criteria used to include certain patients in a follow-up program for potential skin radiation injuries. Peak skin dose resulted the most relevant parameter. The selection criteria refinement (using the cumulative skin dose of the frontal and lateral plane independently and the analysis of the C-arm angulations), together with the introduction of a low-dose protocol in the x-ray system and a training program focused on radiation protection, allowed us to reduce the percentage of patients with doses above trigger levels.

Disclosures: Eliseo Vano—RELATED: Payment for Writing or Reviewing the Manuscript: Private language English reviewer.

REFERENCES

- World Health Organization. *Efficacy and radiation safety in interventional radiology*. Geneva, Switzerland: World Health Organization; 2000
- Valentin J. Avoidance of radiation injuries from medical interventional procedures. *Ann ICRP* 2000;30:7–67
- International Commission on Radiological Protection. Early and late effects of radiation in normal tissues and organs: threshold doses for tissue reactions in a radiation protection context. ICRP Publication 118. *Ann ICRP* 2012;41. In press.
- Stecker MS, Balter S, Towbin RB, et al. Guidelines for patient radiation dose management. *J Vasc Interv Radiol* 2009;20:S263–73
- European Union. Council Directive 97/43 Euratom, on health protection of individuals against the dangers of ionizing radiation in relation to medical exposure, and repealing Directive 84/466 Euratom. Available at: <http://eur-lex.europa.eu/LexUriServ/LexUriServ.do?uri=OJ:L:1997:180:0022:0027:EN:PDF>. Accessed June 9, 2012
- Balter S, Hopewell JW, Miller DL, et al. Fluoroscopically guided interventional procedures: a review of radiation effects on patients' skin and hair. *Radiology* 2010;254:326–41
- Mamourian AC, Young H, Stiefel MF. Cumulative radiation dose in patients admitted with subarachnoid hemorrhage: a prospective study using a self-developing film badge. *AJNR Am J Neuroradiol* 2010;31:1787–90
- Moritake T, Matsumaru Y, Takigawa T, et al. Dose measurement on both patients and operators during neurointerventional procedures using photoluminescence glass dosimeters. *AJNR Am J Neuroradiol* 2008;29:1910–17
- Suzuki S, Furui S, Matsumaru Y, et al. Patient skin dose during neuroembolization by multiple-point measurement using a radiosensitive indicator. *AJNR Am J Neuroradiol* 2008;29:1076–81
- D'Ercole L, Mantovani L, Thyriion FZ, et al. A study on maximum skin dose in cerebral embolization procedures. *AJNR Am J Neuroradiol* 2007;28:503–07
- Moskowitz SI, Davros WJ, Kelly ME, et al. Cumulative radiation dose during hospitalization for aneurysmal subarachnoid hemorrhage. *AJNR Am J Neuroradiol* 2010;31:1377–82
- Miller DL, Balter S, Cole PE, et al. RAD-IR study. Radiation doses in interventional radiology procedures: the RAD-IR study. Part I: overall measures of dose. *J Vasc Interv Radiol* 2003;14:711–27
- D'Ercole L, Thyriion FZ, Bocchiola M, et al. Proposed local diagnostic reference levels in angiography and interventional neuroradiology and a preliminary analysis according to the complexity of the procedures. *Phys Med* 2012;28:61–70
- International Electrotechnical Commission. *Medical electrical equipment. Part 2–43: particular requirements for the safety of x-ray equipment for interventional procedures*. Report IEC 60601-2-43. Geneva, Switzerland: International Electrotechnical Commission; 2010
- Vano E, Geiger B, Schreiner A, et al. Dynamic flat panel detector versus image intensifier in cardiac imaging: dose and image quality. *Phys Med Biol* 2005;50:5731–42
- Martinez LC, Vano E, Gutierrez F, et al. Patient doses from fluoroscopically guided cardiac procedures in pediatrics. *Phys Med Biol* 2007;52:4749–59
- Shope TB. Radiation-induced skin injuries from fluoroscopy. *Radiographics* 1996;16:1195–99
- Vañó E, Arranz L, Sastre JM, et al. Dosimetric and radiation protection considerations based on some cases of patient skin injuries in interventional cardiology. *Br J Radiol* 1998;71:510–16
- Wagner LK, Eifel PJ, Geise RA. Potential biological effects following high x-ray dose interventional procedures. *J Vasc Interv Radiol* 1994;5:71–84
- Koenig TR, Mettler FA, Wagner LK. Skin injuries from fluoroscopically guided procedures: part 2, review of 73 cases and recommendations for minimizing dose delivered to patient. *AJR Am J Roentgenol* 2001;177:13–20
- Mooney RB, McKinstry CS, Kamel HA. Absorbed dose and deterministic effects to patients from interventional neuroradiology. *Br J Radiol* 2000;73:745–51
- Magrassi L, Bongetta D, D'Ercole L, et al. Neuroembolization may expose patients to radiation doses previously linked to tumor induction. *Acta Neurochir* 2012;154:33–41
- Ten JJ, Fernandez JM, Vano E. Automatic management system for dose parameters in interventional radiology and cardiology. *Radiat Prot Dosimetry* 2011;147:325–28
- International Commission on Radiological Protection. Education and training in radiological protection for diagnostic and interventional procedures. ICRP Publication 113. *Ann ICRP* 2009;39:7–68.
- EC Directorate General Environment. Guidelines on education and training in radiation protection for medical exposures. *Radiat Prot* 2000;116. Available at: http://ec.europa.eu/energy/nuclear/radiation_protection/doc/publication/116.pdf. Accessed February 25, 2012
- International Commission on Radiological Protection. Radiological protection in medicine. ICRP Publication 105. *Ann ICRP* 2007;37:1–63.
- Proposal for a Council Directive laying down basic safety standards for protection against the dangers arising from exposure to ionising radiation. Draft presented under Article 31 Euratom Treaty for the opinion of the European Economic and Social Committee. Brussels, Belgium; September 9, 2011. Available at: http://ec.europa.eu/energy/nuclear/radiation_protection/doc/com_2011_0593.pdf. Accessed June 9, 2012